



Skype: chinacaremed

Facebook: China care medical

Wechat: joseph20130708


Whatsapp: +86 139 2648 0041

Email: info@chinacaremedical.com

Website: www.chinacaremedical.com

Add: Room C208, Changsheng Business Building, Daguan South Road No 26, Tianhe, Guangzhou, China

Model CCLQ2508 Rotary Microtome

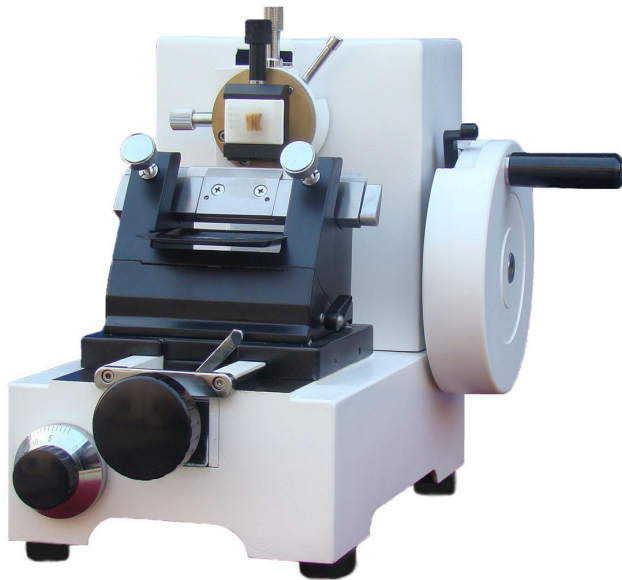

Model CCLQ2508 Rotary Microtome is an improved version on the base of our model CCLQ2508 at the requests of our customers. This model features light feeling, easy and nimble operation, continuous cutting and without the need to move the blade to move or turn the blade holder. This is an ideal instrument for biological histology and clinic pathology.

Main technical performance:

- ◆ The cutter seat can move altering left or right
- ◆ The cutter seat has the function of adjusting angle.
- ◆ The slider of the cutter seat is flexible (using the cross roller guide).
- ◆ Hand wheel has brake function

Main Technical Specifications:

- ◆ Range of Cutting Thickness: 1-40 (μm)
- ◆ Minimum Division of Cutting Thickness Adjustment: 1 (μm)
- ◆ Accuracy of Cutting: $\pm 10\%$



Skype: chinacaremed

Facebook: China care medical

Wechat: joseph20130708

Whatsapp: +86 139 2648 0041

Email: info@chinacaremedical.com

Website: www.chinacaremedical.com

Add: Room C208, Changsheng Business Building, Daguan South Road No 26, Tianhe, Guangzhou, China

- ◆ Maximum Cutting Section: 25×30 (mm)
- ◆ Overall Dimensions: 450×320×320 (mm)

Functions :

- ◆ knob fast trimming, Tumor section, lumber section, cytological section, frozen section

Applications:

- ◆ Its products have been applied to hospitals, universities, research organizations and laboratories involving multiple fields such as biomedicine, life science, agriculture and forestry, health care, disease control and prevention, animal husbandry, and veterinary medicine.